



Mobile Independence Pty Ltd

ABN: 90 531 747 431

1/60 Kremzow Rd
Brendale Qld 4500

Ph: 07 3881 3868
Fax: 07 3881 0525




QUALITY ASSURED COMPANY BY AS/NZS ISO 9001:2008

REGN. Number 3553

Email: sales@mobileindependence.com.au

Web: www.mobileindependence.com.au

General Cushion

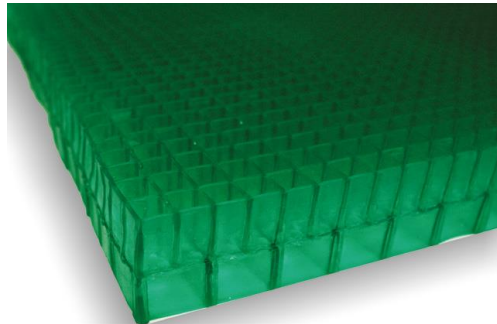

Moderate Pressure Therapy.

Features

- 100% EquaGel™
- 2" tall EquaGel™ in two stages.
- Top layer is thin wall gel and a tight grid for a comfort layer.
- Bottom layer is thicker wall gel and a wide grid for deep posture support.
- Pelvic area of the cushion is contoured to more evenly distribute a user's weight than flat cushions, while retaining 1 1/4" of buckling gel under the pressure points.
- Fluid resistant, vapour permeable, flame resistant cover

Specifications

Sizes	16 X 16
	18 X 16
	18 X 18
	20 X 16
	20 X 18

Warranty: 1.5 year limited