

Pathrider 10 Deluxe

SUPREME PERFORMANCE,
EXHILERATING STYLE,
AND EXCELLENT VALUE.

PATHRIDER 10 DELUXE


The sleek, sporty Pathrider 10 Deluxe scooter delivers smooth performance and all-new features. With a full complement of advanced standard features including, rugged suspension, high-visibility automatic front and rear lighting system, Feather-touch disassembly, LED kerb light, and wraparound easy-drive tiller, the Pathrider 10 Deluxe is a greater value than ever before.



Pathrider 10 Deluxe

Weight capacity:	159 kg (350 lbs.)
Maximum speed ¹ :	Up to 8.5 km/h (5.28 mph)
Ground clearance ² :	3.5 cm (1.38")
Turning radius:	137.16 cm (54")
Overall length ² :	118.75 cm (46.75")
Overall width ² :	56.5 cm (22.24")
Front tyres:	9 cm x 25 cm (3.5" x 10")
Rear tyres:	10 cm x 25 cm (4" x 10")
Suspension:	Front and rear
Range per charge ¹ :	Up to 33.38 km (20.74 miles) with 40 AH batteries
Maximum obstacle climbing ability:	5 cm (2")
Total Weight Without Batteries:	84 kg (185 lbs.)
Heaviest piece when disassembled:	29.48 kg (65 lbs.)
Standard seating:	Type: CRS Junior seat Weight: 19 kg (42 lbs.) Material: Black Vinyl Dimensions: width – 44.5 cm (17.5") depth – 43 cm (17")
Drive System:	Rear-wheel drive, sealed transaxle, 24-volt DC motor
Dual Braking System:	Electronic, regenerative, and electromechanical
Battery Requirements ³ :	Two 12-volt, deep cycle Size: 35 AH or optional 40 AH Weight ⁴ : 35 AH: 11 kg (24.5 lbs.) each; 40 AH: 15 kg (33 lbs.) each
Battery Charger	5-amp, off-board

Colours:

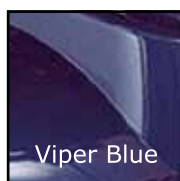

¹Varies with user weight, terrain type, battery amp-hour (AH), battery charge, battery condition and tyre condition.

²Due to manufacturing tolerances and continual product improvement, this specification can be subject to a variance of (+ or - 3%).

³AGM or Gel-Cell type required.

⁴Battery weight may vary based on manufacturer.

NOTE: All specifications subject to change without notice.

Features

- Front and rear suspension
- High-visibility automatic front and rear lighting system
- Exclusive one-hand Feather-touch disassembly provides the ultimate quick transport system
- Long-lasting, high-intensity LED headlight, tail light and indicators
- Pride's exclusive stylish, light weight, non-marking, low profile wheels
- Auto-connecting front-to-rear harness
- Easy-drive tiller with wraparound handles
- New back-lit battery meter for increased accuracy and visibility in any setting
- Includes pinchless hinges & eight-position seat rotation
- Easy-grip tiller adjustment knob
- Easy access tiller-mounted charger port for standard dual voltage off-board charger
- Rear view mirror standard

Options

- Forearm crutch holder
- Cane/crutch holder
- Oxygen tank holder
- Cup holder
- Rear basket
- Safety flag
- Walker holder
- Weather cover
- Saddlebag

The information contained herein is correct at the time of publication; due to Pride's commitment to constant improvement and development, we reserve the right to alter specifications without prior notice. Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressures, ambient temperature, and the gradient and surface of the road or pavement.

